

Bruno Elite Curve

Custom Curved Rail Stairlifts

*Turning
Houses...
Back Into
Homes!*



 **BRUNO**[®]
for your independence

 **MADE IN USA**

Superior Styling



Your home's lovely curved staircase was once a highlight. Now, it's become a challenge. Bruno's Custom Curved Rail Stairlifts restore the function of your staircase for everyone, in beautiful, graceful style.

Refinement

This is the benchmark for stairlift design and operation. The Elite Curve, Model CRE-2110, encompasses all the technological advances and aesthetic considerations in one magnificent stairlift. Using two 12V batteries which are powered by a battery charger that plugs into any household outlet, the Elite Curve provides uninterrupted access up and down the stairway, even if your home experiences a power outage.

Convenient Installation

The Elite Curve is delivered fully equipped and ready for installation. The unit anchors to the steps, not the wall, and the stairlift rail can be placed as close as 7" from the wall. When the chair's seat, footrest and armrests are folded, there is plenty of room for family and guests to use the staircase unimpeded.

Bruno's wireless call/send system allows for faster and cleaner installations. Your custom designed rail has been assembled and tested extensively at Bruno's manufacturing facilities, ensuring a perfect fit in your home.



When the Elite Curve is in the folded position, with seat, footrest and armrests raised, there is plenty of room for family and guests to use the stairway.



The simple, unique design of the rail hides all mechanical components, and installation requires no modification to your staircase.

Elite Curve


Quality Craftsmanship

Smooth & Stable

The Elite Curve delivers an effortless ride with subtle power and remarkable stability over any stairway configuration. Its strong, resilient drive system provides the safest operation and longest lasting durability for years to come, even when carrying a maximum load of 400 lb. And, the soft-start, soft-stop, solid state controller features a secure, comfortable ride.

Designed for You!

Every aspect of the Elite Curve is a work of art. The seat, with upholstery selections, and carriage, incorporates everything Bruno has learned in the creation of the world's most comfortable, reliable stairlifts. The elegant, curved rail is produced by experienced craftsmen from exact measurements of your staircase. No two stairlifts are alike. Only the quality is consistent.

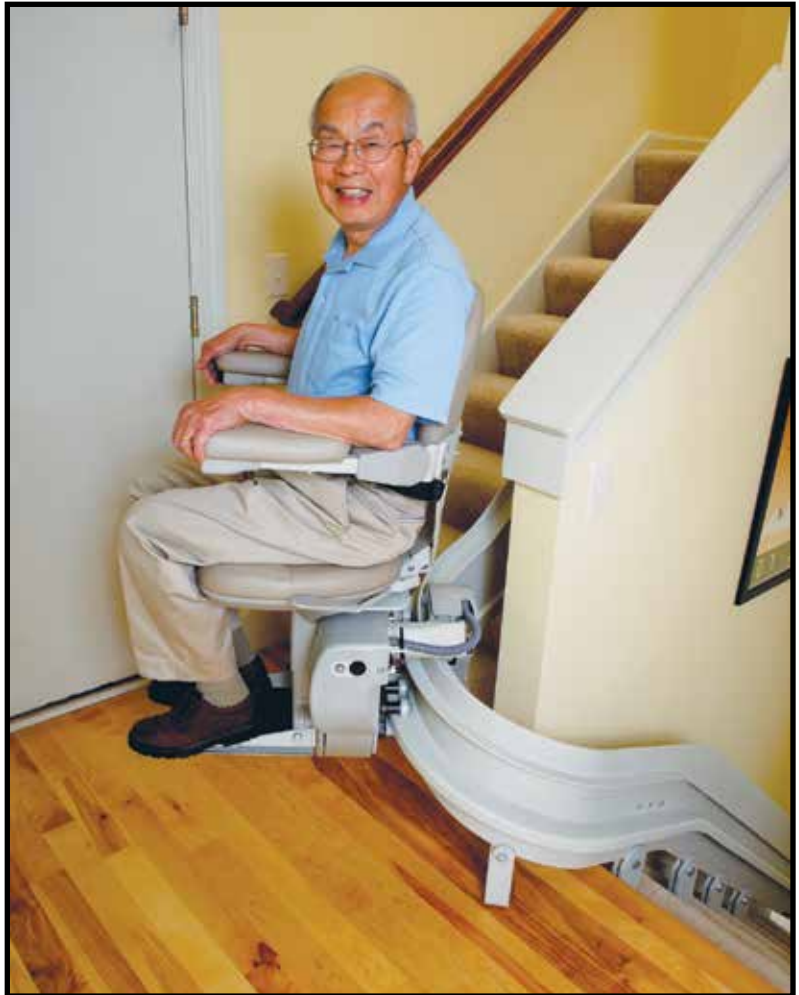

Elite Curve Features



The offset swivel seat makes the entry and exit position at the top landing safer and easier than ever!



The footrest and carriage safety sensors safely stop you, even if the slightest obstacle is in the way while you ride!



RAIL OPTION:
The CRE-2110 can be equipped with an optional top or bottom park to facilitate an exit and remove the unit from the stairs.



The two wireless call/send controls make installation simple and clean with no wires running along the wall.

Elite Curve

Elite Curve CRE-2110 Specifications

Exceptional Value

- ◆ Locking swivel seat
- ◆ Two remote, wireless call/send controls
- ◆ Generous seat size – 17.75" (451 mm) adjustable to 21.75" (552 mm) between the arms
- ◆ Footrest and carriage safety sensor and retractable seat belt
- ◆ Charging-station indicator lights at the top and bottom of the stairway
- ◆ Onboard audio and visual diagnostics make service easy and precise
- ◆ Meets all ANSI - A18.1 requirements; CSA/C-US listed
- ◆ Bruno's Gold Warranty – Five years on major components and two years on parts

Convenient Installation

- ◆ Rail typically installs to within 7" (178 mm) of the wall - varies with each application
- ◆ Installations can be done on either the inside or outside of most staircases.

Comfort

- ◆ Adjustable seat height for easier entry/exit at the top and bottom of the stairway
- ◆ Adjustable footrest height for added comfort
- ◆ Flip-up arms for easy wheelchair transfers
- ◆ Adjustable armrest width
- ◆ Contoured seat for maximum comfort

Battery-Powered Technology

- ◆ Ensures access even during power outages
- ◆ Powered by two 12V batteries; continuously charged by a 2A battery charger that plugs into any 110V household outlet
- ◆ Requires no special wiring

Smooth, Powerful Drive System

- ◆ Self-locking worm gear; rack-and-pinion drive
- ◆ Maximum rider weight of 400 lb (181 kg)
- ◆ Soft-start, soft-stop for the user's maximum comfort

Ask us about our other...

Home accessibility solutions and automotive products



OUTDOOR ELITE Model SRE-2010E
The only 400 lb/181.8 kg capacity outdoor stairlift... an INDUSTRY EXCLUSIVE!

VERTICAL PLATFORM LIFT
Model VPL-3100

VEHICLE LIFTS

VALET® SIGNATURE SEATING

Founded over 25 years ago, Bruno Independent Living Aids, Inc., is a world-wide manufacturer of accessibility products designed to enhance the lives of those challenged by limited mobility. In pursuing this goal, Bruno is proud to be the only U.S. manufacturer of indoor/outdoor stairlifts, interior/exterior scooter, powerchair and wheelchair vehicle lifts and Valet® Signature Seating products to be ISO 9001 Certified. This is one of the most stringent, widely-recognized quality standards in the world and your assurance that Bruno is committed to providing you with the highest-quality, most reliable products in the industry. Ask your certified Bruno dealer for details on all Bruno independence products.

Represented by:



BRUNO INDEPENDENT LIVING AIDS, INC.®

1780 Executive Drive, Oconomowoc, WI 53066

www.bruno.com



All illustrations and specifications in this brochure are based on the latest product information available at the time of publication. Bruno Independent Living Aids, Inc. reserves the right to make changes at any time without notice. © Bruno Independent Living Aids, Inc. 2013. 2013-9-4_Elite Curve CRE-2110.indd

Part #: MKT-CRE-2110